




RELY ON PRECISION – **PANASIL®**

A-SILICONE FOR PRECISION IMPRESSION TAKING

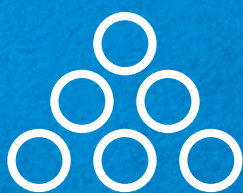

KETTENBACHDENTAL


Simply intelligent

PRECISION WITHOUT COMPROMISES

Precision impression material (A-silicone)

In the **Panasil®** product family you find perfectly coordinated, innovative materials for all impression techniques and indications.

DOUBLE-MIX IMPRESSIONS

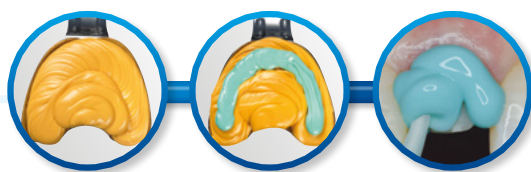

Panasil® tray Soft Heavy

PLEASANTLY HEAVY-BODIED



Panasil® initial contact Light




Panasil® contact two in one Light

LOW VISCOSITY

INGENIOUSLY HYDROPHILIC

FAST SETTING TIME

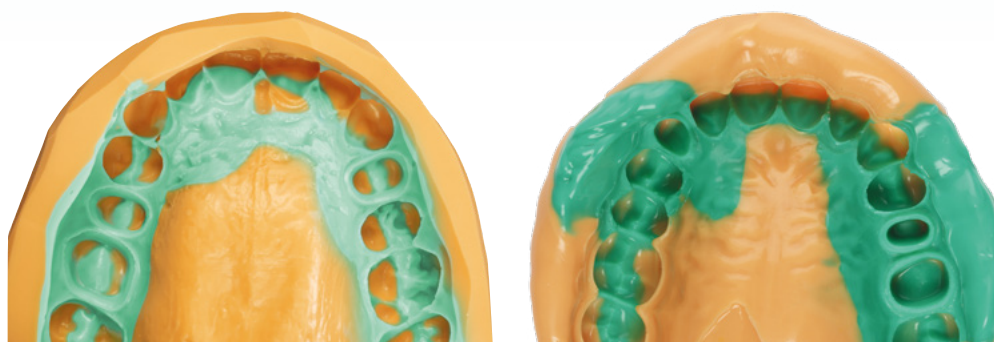

FOR EVERY TECHNIQUE AND EVERY INDICATION

TWO-STEP IMPRESSIONS



THE FAVOURITE



Panasil® binetics Putty Fast

EASY TO TRIM



Panasil® Putty Fast

**PRECISE AND
ECONOMICAL**



Panasil® tray Fast Heavy

FAST SETTING TIME



Panasil® initial contact X-Light




Panasil® contact plus X-Light

VERY LOW VISCOSITY

INGENIOUSLY HYDROPHILIC

FAST SETTING TIME

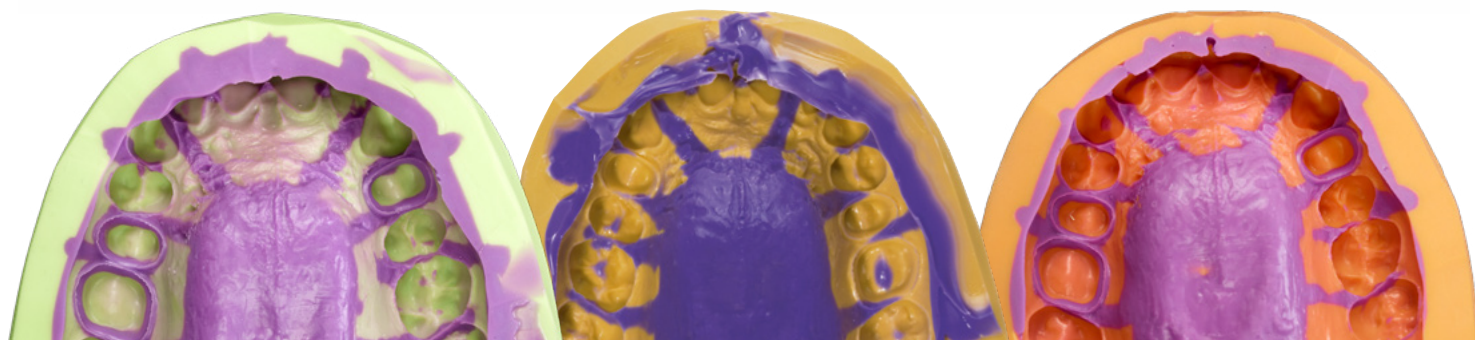

SANDWICH IMPRESSIONS



Panasil®
binetics Putty Soft


Panasil® Putty Soft

EASY TO KNEAD, STABLE



Panasil® initial contact Light

**LOW VISCOSITY,
INGENIOUSLY HYDROPHILIC**

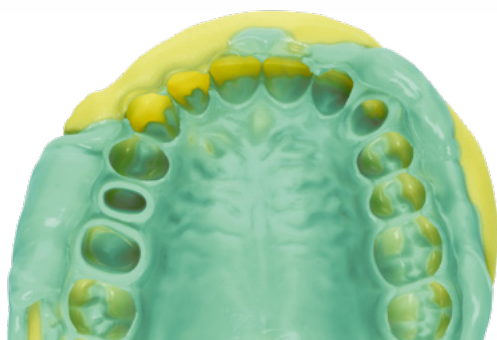

MONOPHASE IMPRESSIONS



Panasil® monophase Medium

HIGH FINAL HARDNESS



Panasil® monophase Medium

HIGH FINAL HARDNESS




OUTSTANDINGLY HYDROPHILIC FROM THE VERY FIRST SECOND

As the contact angle measurement clearly confirms

(parameter for the hydrophilicity of a material)

Panasil® initial contact reaches the lowest contact angle after just three seconds compared to the competitor products mentioned.

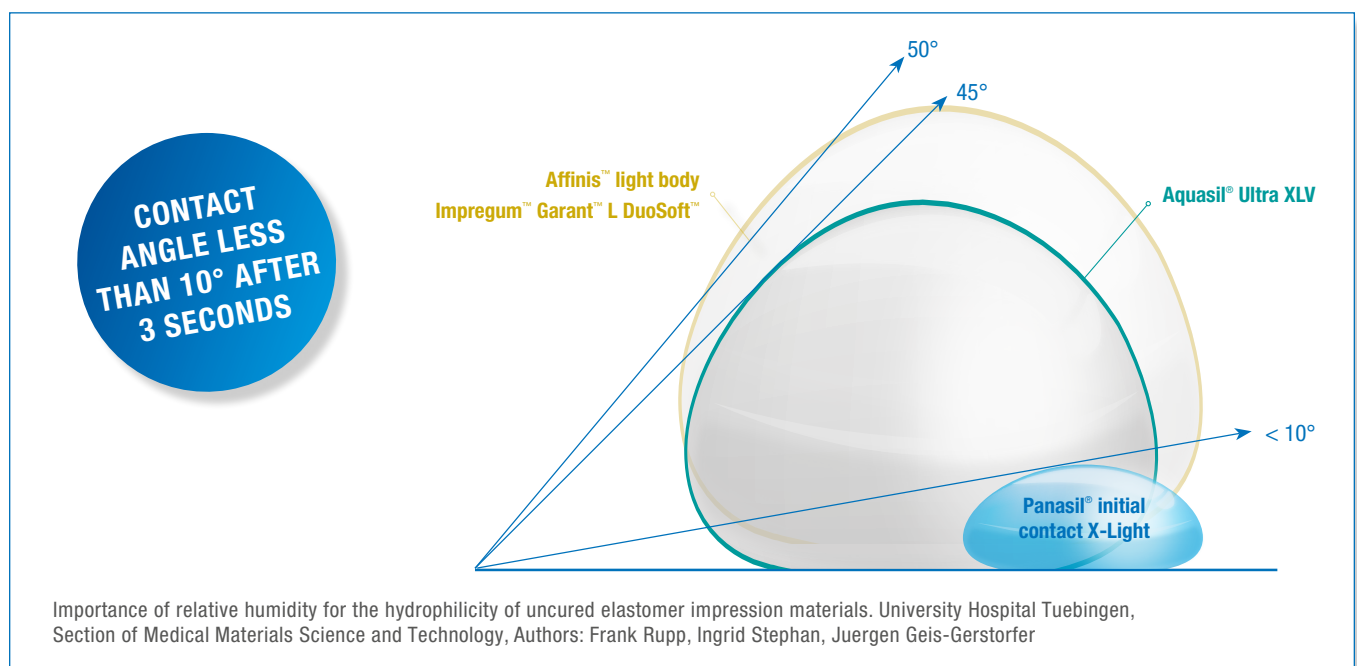

The flatter the drop, the more pronounced the hydrophilicity of the material.

EXCELLENT HYDROPHILICITY



- **Panasil® immediately overcomes** any **residual moisture** on direct contact with the tooth surface and ensures a precise result.



- **Panasil® ensures** sharply defined preparation margins even in **extreme situations**, e.g. in the case of difficult hemostasis.

PANASIL® – THE RANGE

PUTTY	Panasil® Putty		Panasil® Putty Soft	Panasil® Putty Fast
	Intro pack	–	REF 11120	REF 11140
		200 mL base paste, 200 mL catalyst paste, 2 dispensing scoops		
	Normal pack	REF 11101	REF 11121	REF 11141
		450 mL base paste, 450 mL catalyst paste, 2 dispensing scoops		
	Economy pack	REF 11103	REF 11123	REF 11143
		4 x 450 mL base paste, 4 x 450 mL catalyst paste, 2 dispensing scoops		
	Panasil® binetics Putty Soft		Panasil® binetics Putty Fast	
	Refill pack	REF 14703		REF 14701
		2 x 380 mL cartridges		

HEAVY	Panasil® tray Soft Heavy		Panasil® tray Fast Heavy	Panasil® tray Soft Heavy Fast
	Refill pack	REF 14707	REF 14705	–
		2 x 380 mL cartridges		
	Normal pack	REF 13541	REF 13551	REF 13561
		2 x 50 mL cartridges, 6 white mixing tips		

MEDIUM	Panasil® monophasic Medium			
	Refill pack		REF 14709	
		2 x 380 mL cartridges		
	Normal pack		REF 13501	
		2 x 50 mL cartridges, 6 white mixing tips		

LIGHT	Panasil® initial contact X-Light/ X-Light Fast		Panasil® initial contact Light/ Light Fast	Panasil® initial contact Regular
	Normal pack	REF 13401/ REF 13461	REF 13411 / REF 13471	REF 13431
		2 x 50 mL cartridges, 8 grey mixing tips		
	Bonus pack	REF 28300 / –	REF 28310 / –	–
		10 x 50 mL cartridges, 40 grey mixing tips		
	Panasil® contact plus X-Light		Panasil® contact two in one Light	
	Normal pack	REF 11892		REF 11781
		2 x 50 mL cartridges, 8 grey mixing tips		
			2 x 50 mL cartridges, 6 white mixing tips	



1800 651 788 www.pharmacodental.com.au

ALWAYS READ THE LABEL AND FOLLOW THE DIRECTIONS FOR USE.

Some medical devices are for healthcare professional use only.
Some products or variants may not be registered or available in Australia.

Pharmaco (Australia) Ltd, Chatswood NSW 2067 PHAA443