




SOOO ADHERENT IS **MUCOPREN®** SOFT

SOFT RELINING

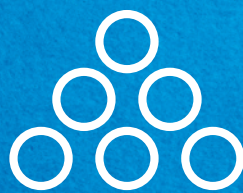


KETTENBACHDENTAL
Simply intelligent

MUCOPREN® SOFT

Prosthetic relining

Mucopren® Soft is a soft long-term relining material based on a silicone for **direct and indirect application**.

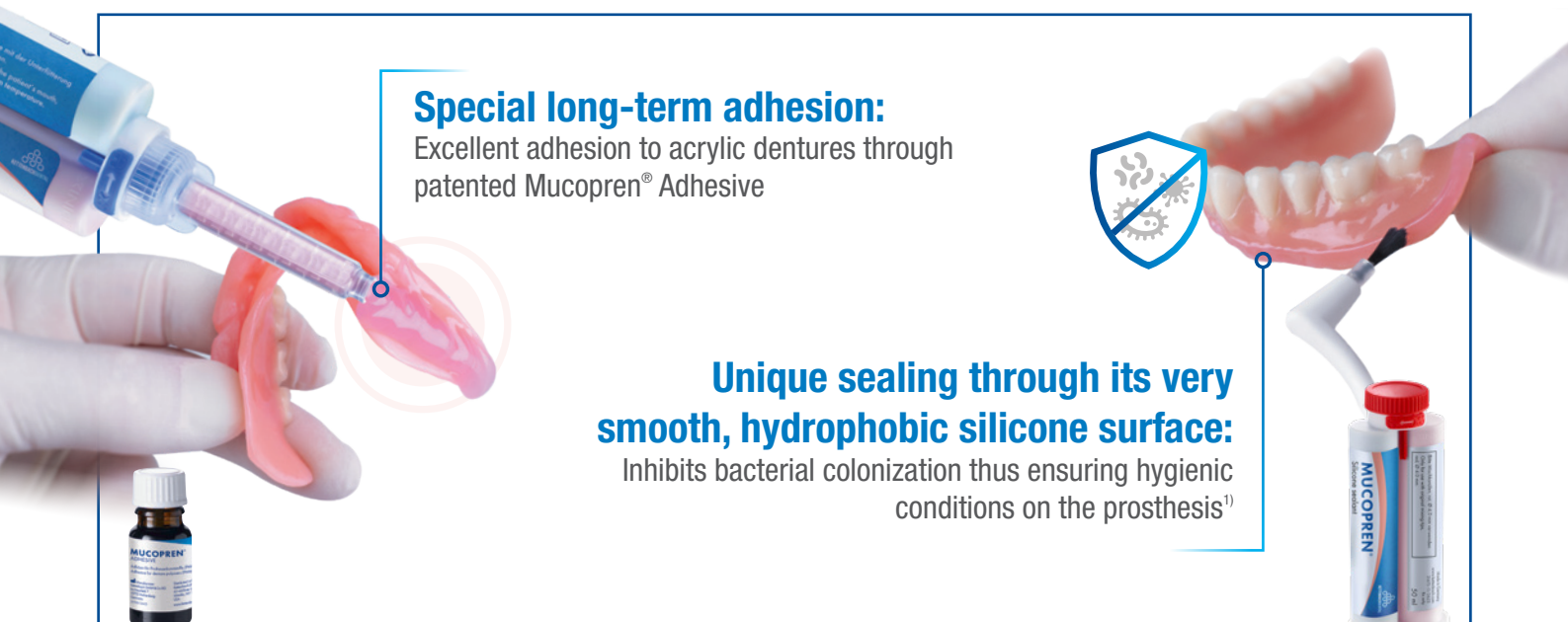

Special long-term adhesion:
Excellent adhesion to acrylic dentures through patented Mucopren® Adhesive

Unique sealing through its very smooth, hydrophobic silicone surface:
Inhibits bacterial colonization thus ensuring hygienic conditions on the prosthesis¹⁾



Application video

Simple treatment & fast processing

Faster setting times and shorter waiting times ensure a fast and smooth process in the practice and laboratory.

Supports the healing process of implants

Also suitable for use in implantology for trouble-free osseointegration.



High wearing comfort

The definitive resilient reline is completely odorless and tasteless and makes wearing the prosthesis as pleasant as possible.



¹⁾ Result from in-vitro studies in which it was shown that the growth of *Candida albicans* was prevented thanks to the smooth silicone surface

HOLDS TO ITS PROMISES



FURTHER INDICATIONS



- For the creation of **permanently soft relinings**
- For relieving pressure on the surgical site after insertion of **implants**
- For balancing mucous membrane areas with different resilience
- For relieving sharp bony ridges
- For relining **interim prostheses**
- As **mucous membrane conditioner** in case of prosthetic stomatitis
- For creating and relining **obturators**

Benefits of using the new adhesive

» ONLY A **SINGLE APPLICATION** NECESSARY



3 MIN

CHAIRSIDE

The processing time is reduced to only **3 minutes** following mouth removal. No water bath necessary.



3 MIN

LABSIDE

The processing time is reduced to only **3 minutes** in the water bath or pressure pot.

MUCOPREN® SOFT – THE RANGE



Mucopren® Soft

Mucopren® Soft Basic Set	REF 28105	50 mL Mucopren® Soft, 50 mL silicone sealant, 10 mL adhesive, 6 blue mixing tips, 16 red mixing tips, 1 brush holder, 20 single-use brushes, 1 steel bur, accessories	
Mucopren® Soft Normal pack	REF 15687	2 x 50 mL Mucopren® Soft, 6 blue mixing tips	
Mucopren® Silicone sealant	REF 15686	50 mL Mucopren® Silicone sealant, 8 red mixing tips	
Mucopren® Adhesive	REF 14203	10 mL Mucopren® Adhesive	



Technical data

Product name	Mixing ratio	Working time at 73 °F ≤	Intraoral setting time at 95 °F ≥	Water bath	End of setting ≥	Shore hardness (approx.)
Mucopren® Soft – direct relining	1:1	2:15 min.	3:15 min.	Not applicable	5:30 min.*	A 28
Mucopren® Soft – indirect relining	1:1	2:15 min.	Not applicable	3:00 min.	> 3 min.	A 28
Mucopren® Soft silicone sealant	1:1	1:30 min.	Not applicable	Not applicable	6:30 min. at 23°C	A 33

*Total setting time (removal time from mouth) from start of mixing



1800 651 788 www.pharmacodental.com.au

ALWAYS READ THE LABEL AND FOLLOW THE DIRECTIONS FOR USE.

Some medical devices are for healthcare professional use only.
Some products or variants may not be registered or available in Australia.